


Pharmacare Health Specialists

2740 Regency Park Wichita, KS 67226

Phone 316-681-2181 Fax 316-681-0277

Toll Free 800-253-9411

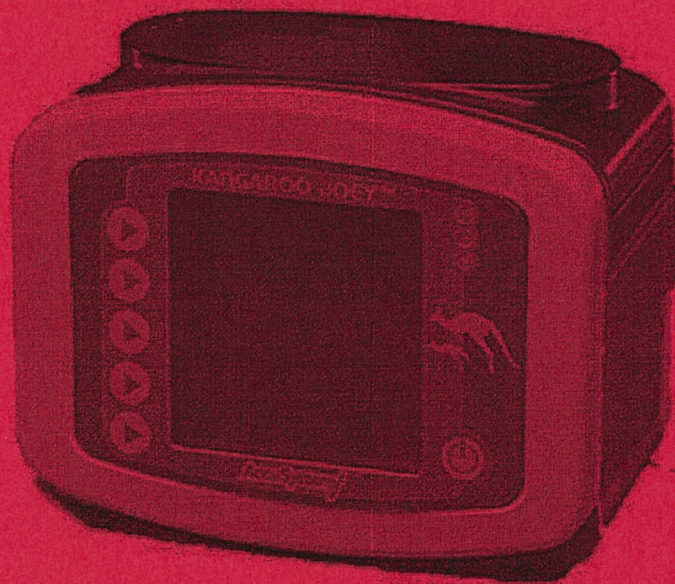

Kangaroo Joey Instructions

KENDALL

KANGAROO JOEY™


Quick Reference Guide

Additional Operations


- For additional operations or information, see the KANGAROO JOEY enteral feeding pump Operation and Service manual.
- See back for error informational conditions.




General Startup

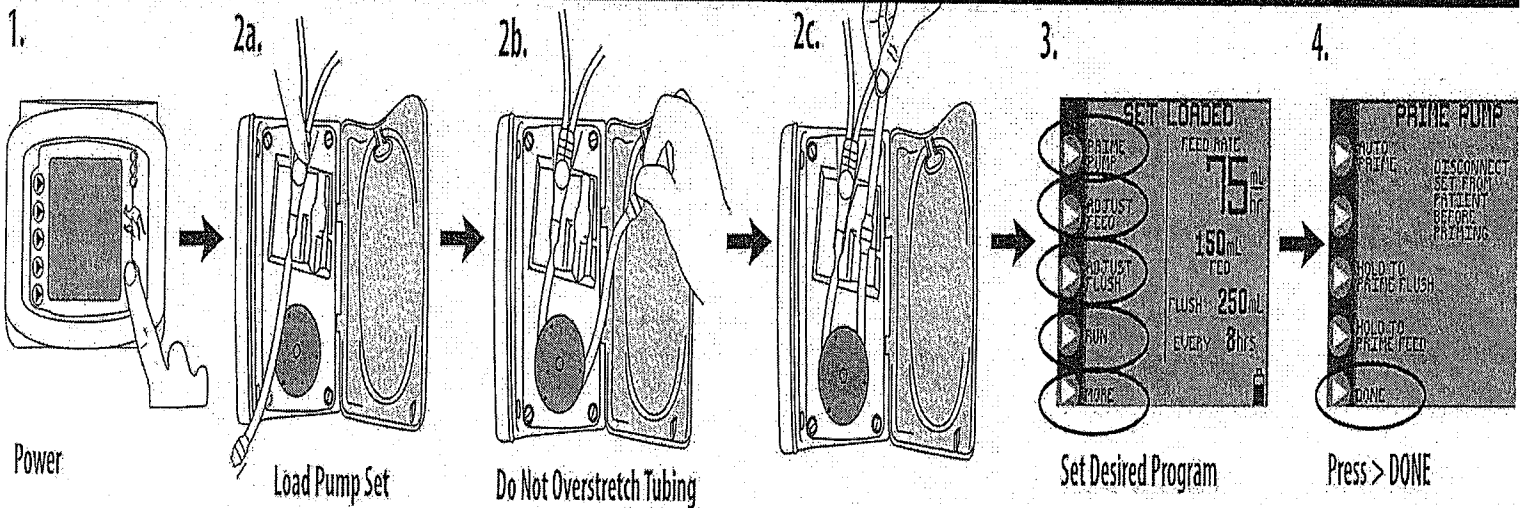

Priming Options

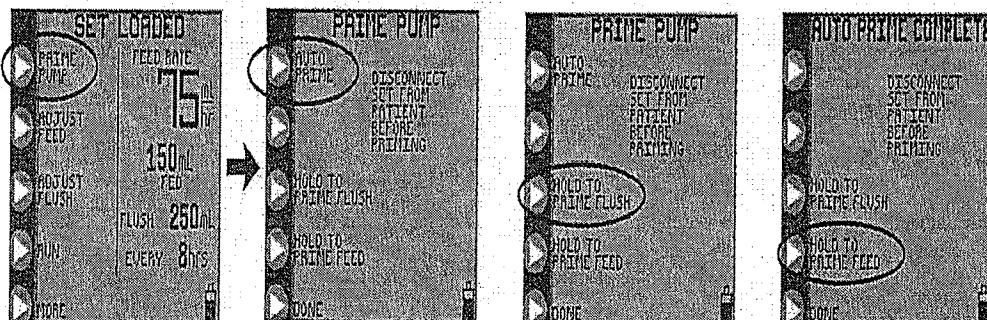

Program Continuous Mode

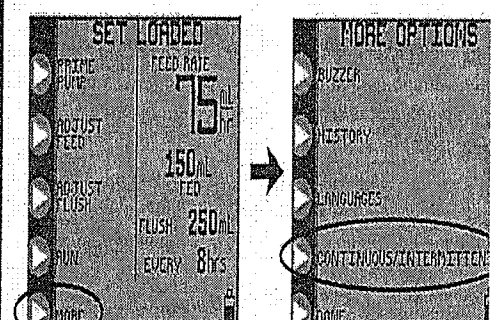

™ Trademark of Tyco Healthcare Group LP or its affiliate.

Patents Pending

TYCO HEALTHCARE GROUP LP • 15 HAMPSHIRE STREET • MANSFIELD, MA 02048 • © 2008 TYCO HEALTHCARE GROUP LP
ALL RIGHTS RESERVED • 080625 • PRODUCT INFORMATION IN U.S. 1-800-962-9888 • www.tycohealthcare.com

EC REP TYCO HEALTHCARE U.K. LTD. • GOSPORT PO13 0AS, U.K.



0123

Program Intermittent Mode

1. **SET LOADED**
PRIME PUMP
FEED RATE 75 mL/hr
ADJUST FEED 150 mL FED
ADJUST FLUSH 250 mL
RUN EVERY 8 hrs
MORE

1a. **MORE OPTIONS**
BUZZER
HISTORY
LANGUAGES
CONTINUOUS/INTERMITTENT
DONE

2. **SET LOADED**
PRIME PUMP
FEED RATE 400 mL/hr
ADJUST FEED 375 mL
ADJUST FLUSH 375 mL
RUN EVERY 6 hrs
FLUSH 250 mL
EVERY 8 hrs
MORE

Program Flushing

SET LOADED
PRIME PUMP
FEED RATE 75 mL/hr
ADJUST FEED 150 mL FED
ADJUST FLUSH 250 mL
RUN EVERY 8 hrs
MORE

Clear Volume Delivered

RUNNING
FEED RATE 75 mL/hr
ADJUST FEED 150 mL FED
ADJUST FLUSH 250 mL
RUN EVERY 8 hrs
HOLD

HOLDING
FEED RATE 75 mL/hr
ADJUST FEED 150 mL FED
ADJUST FLUSH 250 mL
RESUME IN MIN
RUN
HISTORY

Operations

Review Program History

SET LOADED
PRIME PUMP
FEED RATE 75 mL/hr
ADJUST FEED 150 mL FED
ADJUST FLUSH 250 mL
RUN EVERY 8 hrs
MORE

MORE OPTIONS
BUZZER
HISTORY
LANGUAGES
CONTINUOUS/INTERMITTENT
DONE

Clear All Settings

PRIOR SETTINGS
KEEP SETTINGS
CLEAR SETTINGS

At Start Up

Power Off

HOLD BUTTON TO POWER DOWN
2

Error Informational Conditions

HOLD ERROR
CONTINUE
POWER DOWN

Pump HOLDING 10 Minutes

FLOW ERROR
CLOG IN LINE DOWNSTREAM OF PUMP
CONTINUE
STOP ALARM

Occlusion

FEED ERROR
- BAG EMPTY
- CLOG IN LINE
- VALVE NOT LOADED
CONTINUE
STOP ALARM

Empty Feed Bag Occlusion

FLUSH ERROR
- BAG EMPTY
- CLOG IN LINE
- VALVE NOT LOADED
CONTINUE
STOP ALARM

Empty Flush Bag Occlusion

RUNNING
FEED RATE 400 mL/hr
ADJUST FEED 375 mL
ADJUST FLUSH 375 mL
RUN EVERY 6 hrs
HOLD

Pump Set In Use > 24 Hours

PUMP SET DISLOADED
POWER DOWN AND RE-INSTALL SET
POWER DOWN

Pump Set Not Properly Loaded

BATTERY LOW
PLUG IN POWER CORD
CONTINUE
POWER DOWN


Low Battery

ROTOR ERROR
CONTINUE
STOP ALARM

Pump Set Tubing Not Properly Loaded On Rotor, Or Unusual Rotor Operating Condition

SYSTEM ERROR
PLEASE CALL TECH SUPPORT
50
POWER DOWN

General Error. See Manual

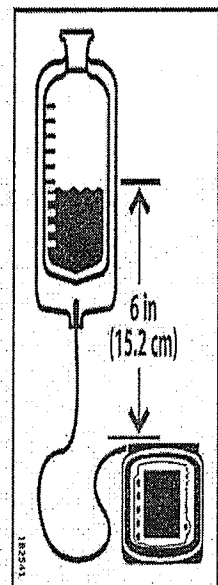

Correct Fluid Level